

1

foin

TP 075 – BANANA MULTIPLICATION MEDIUM (W/Vitamins, Surcose, Tryptone, Agar & IAA)

INTENDED USE

Supports or facilitate plant growth and/or shoot proliferation in two or more plant tissue cultures (both in monocotyledons and dicotyledons).

COMPOSITION

Ingredients	Mg / Ltr
Copper sulphate (Penta hydrate)	0.025
Cobalt chloride (Hexa hydrate)	0.025
Boric acid	6.200
Potassium iodide	0.830
Zinc sulphate (Hepta hydrate)	8.600
Sodium molybdate	0.250
Na2EDTA	37.260
Ferrous sulphate (Hepta hydrate)	27.800
Manganese sulphate (Mono hydrate)	16.900
Sodium phosphate (Mono basic)	221.710
Potassium phosphate (Mono basic)	170.000
Magnesium sulphate	180.690
Calcium chloride (Dihydrate)	440.000
Ammonium nitrate	1650.000
Potassium nitrate	1900.000
Myo-inositol	100.000
Thiamine hydrochloride	0.400
Indole-3-acetic acid	1.000
Casein hydrolysate	10.000
Sucrose	30000.00
Agar	8000.00

Formula weight: 42.77 gms/ltr

: 5.8 ± 0.2

QUALITY CONTROL SPECIFICATIONS


Appearance of Powder	
Appearance of prepared medium	
pH (at 25°C)	

: White to light tan with homogenous mixture of free flowing powder.

: Colorless to slight yellow solution, clear, complete.

INSTRUCTION FOR USE

Dissolve 42.77 gms of dehydrated medium in 600 ml of distilled or deionized water at room temperature (15-30°C). Rinse media vial with small quantity of distilled water to remove traces of powder. Add the desired heat stable supplements prior to autoclaving. Continue stirring until the powder has dissolved. Sometimes media does not dissolve

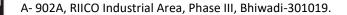


PRODUCT DATA SHEET



completely unless the pH is reduced. For these, lower the pH to about 3.0 to facilitate dissolution of media. The pH of medium is adjusted by using 1N HCL/ 1N NaOH/ 1N KOH. Make up the final volume to 1000ml with distilled water. Mix gently, heat and rotate between intervals until the solution becomes clear. Do not boil, reheat and allow to cool below 50°C during dispensing. Dispense the medium into suitable containers, plug or cap, then autoclave at 15psi (121°C) for 15 minutes, using a slow exhaust cycle. Higher temperatures and/or longer times are not recommended. Cool the autoclaved culture vessels containing medium to 45-50°C and aseptically add desired sterile heat-labile substrate. **Note**: Media should be prepared according to formula mentioned on the label however, it is recommended to use an entire container at once. Heat-labile substrates should be added, after autoclaving.

STORAGE

Dehydrated plant tissue culture media is hygroscopic and should be protected from sunlight and moisture. Store the prepared medium at 2–8°C away from direct light. Medium should be used before the expiry date. The entire volume of each bottle should be used immediately after opening or else the unused portion should be stored in desiccators and refrigerated at 2–8°C.

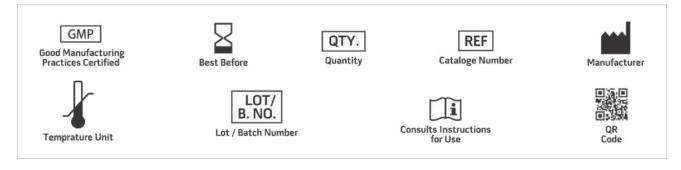

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices. *For Lab Use Only Revision: 28 Feb., 2022

