

TS 021 – HEAMOGLOBIN POWDER

INTENDED USE

For isolation of *Neisseria* spp.

COMPOSITION

Ingredients	Concentration
Haemoglobin powder	50G/ 100G

Applied in Medium

TM 064 - CHOCOLATE AGAR BASE

TM 116 - GC AGAR BASE

TMV 116 - GC AGAR BASE (VEG.)

TM 879 - TELLURITE BLOOD AGAR BASE

TM 933 - THAYER MARTIN MEDIUM BASE

TM 1105 - TRANSGROW MEDIUM BASE


TM 2264 - CHOCOLATE NO. 2 AGAR BASE

INSTRUCTION FOR USE

A specially prepared powder whose 2% w/v aqueous solution is autoclavable. The aqueous solution is chocolate brown, opaque and contains flocculent dispersible precipitate. It is used for 500 ml medium preparation of TM 116 - GC AGAR BASE - 5 gms / TMV 116 - GC AGAR BASE (VEG.) - 5 gms / TM 933 - THAYER MARTIN MEDIUM BASE – 5 gms / TM 2264 - CHOCOLATE NO. 2 AGAR BASE – 5 gms / TM 879 - TELLURITE BLOOD AGAR BASE – 10 gms / TM 064 - CHOCOLATE AGAR BASE - 10 gms / TM 1105 - TRANSGROW MEDIUM BASE - 2 gms.

STORAGE

Vials should be stored in sealed container at 2-8°C.

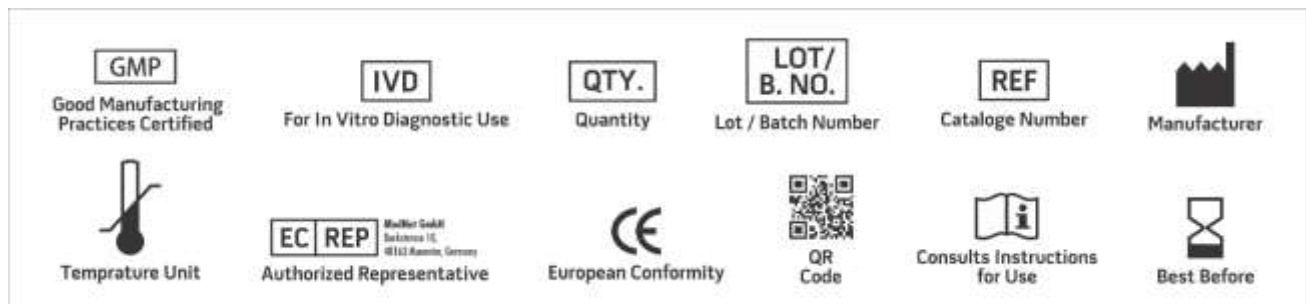

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 2 March., 2022