

TS 089 – STAPH STREPTO SUPPLEMENT

INTENDED USE

For selective isolation of Staphylococci and Streptococci.

COMPOSITION

Ingredients	Concentration
Nalidixic acid	1mg
Colistin sulphate	1mg

(per vial sufficient for 500 ml medium)

Applied in Medium

TM 071 - COLUMBIA BLOOD AGAR BASE

TMV 071 - COLUMBIA BLOOD AGAR BASE (VEG.)


TM 2044 - COLUMBIA BLOOD AGAR BASE W/ 1% AGAR

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial aseptically with 2 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 071 - COLUMBIA BLOOD AGAR BASE / TMV 071 - COLUMBIA BLOOD AGAR BASE (VEG.) along with 25 ml of sterile defibrinated horse blood, TM 2044 - COLUMBIA BLOOD AGAR BASE W/ 1% AGAR. Mix gently and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

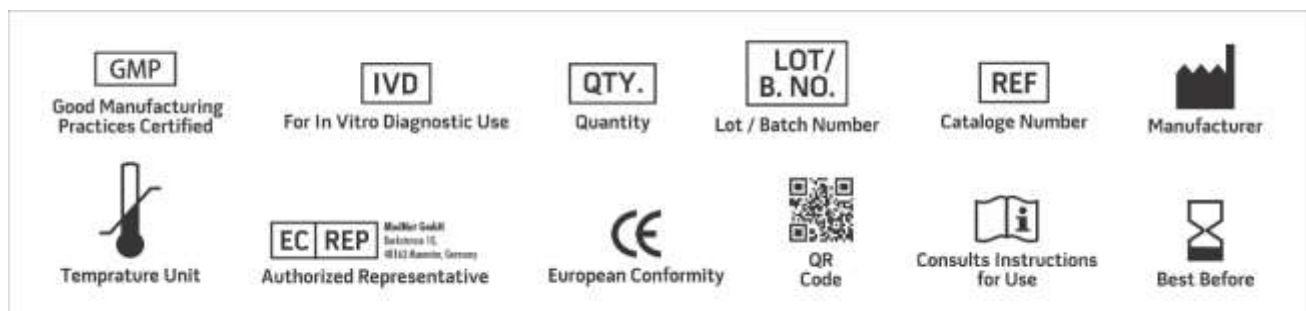

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 3 March., 2022