

TS 108 – SODIUM BICARBONATE SOLUTION (20ml/vl)

INTENDED USE

For selective isolation of *Enterococcus* species.

COMPOSITION

Ingredients	Concentration
Sodium bicarbonate (10%w / v)	20ml

(per vial sufficient for 980 ml medium)

Applied in Medium


TM 1190 - GTC AGAR BASE

INSTRUCTION FOR USE

Warm up the refrigerated solution to room temperature. Add content of 1 vial aseptically with 980 ml sterile, molten, cooled (45-50°C) TM 1190 - GTC Agar Base M922 along with 2 vials of GTC Supplement. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

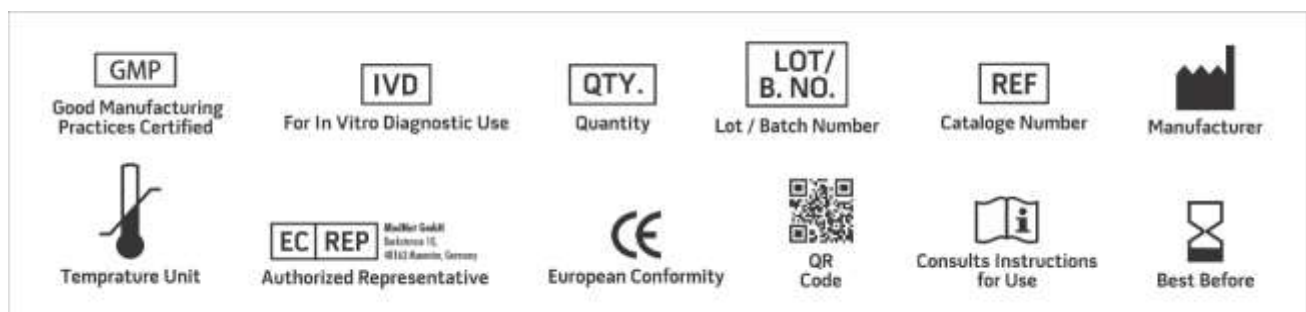

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

*For Lab Use Only
Revision: 3 March., 2022