

TS 110 – KIMMIG SELECTIVE SUPPLEMENT

INTENDED USE

For selective isolation and cultivation of fungi.

COMPOSITION

Ingredients	Concentration
Colistin sulphate	40mg
Novobiocin	50mg

(per vial sufficient for 500 ml medium)

Applied in Medium


TM 1211 - KIMMIG FUNGI AGAR BASE

INSTRUCTION FOR USE

Rehydrate the content of one vial aseptically with 2 ml of sterile distilled water. Mix well and aseptically add it to 500 ml of sterile, molten, cooled (45-50°C) TM 1211 - KIMMIG FUNGI AGAR BASE. Mix well and aseptically pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

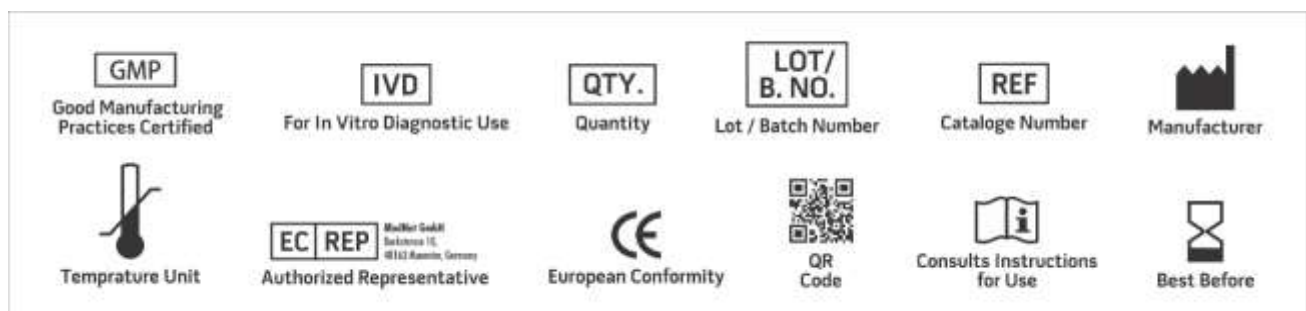

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 2 March., 2022