

TS 161 – LACHICA SUPPLEMENT

INTENDED USE

for selective isolation of *Aeromonas hydrophilla*.

COMPOSITION

Ingredients	Concentration
Ampicillin	10mg

(Per vial sufficient for 1000 ml medium)

Applied in Medium


TM 1369 - LACHICA'S MEDIUM BASE

INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 2 ml of 50% v/v alcohol. Mix well and aseptically add to 1000 ml sterile, molten TM 1369 - LACHICA'S MEDIUM BASE. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

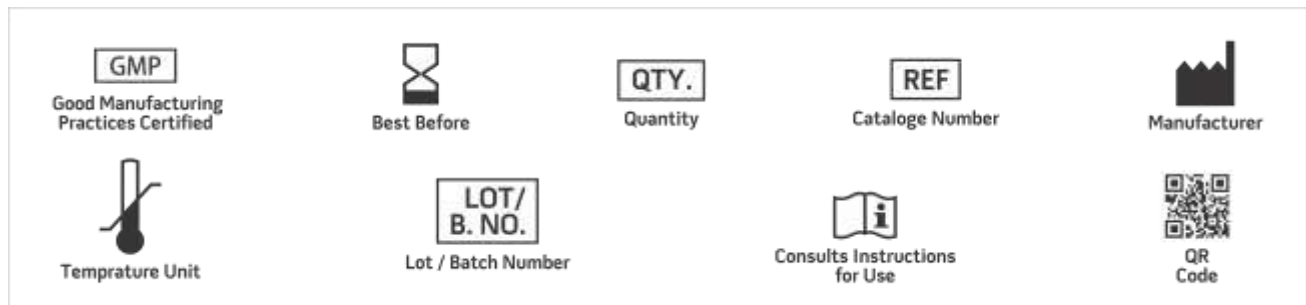

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 7 March., 2022