

## TS 218 – VANCOMYCIN SUPPLEMENT

### INTENDED USE

For isolation of *Cronobacter sakazakii*.

### COMPOSITION

Ingredients	Concentration
Vancomycin	3mg

(per vial sufficient for 500 ml medium)

### Applied in Medium

TM 1642 - VRE AGAR BASE (VANCOMYCIN RESISTANT ENTEROCOCCI AGAR)


TM 2047 - CRONOBACTER SELECTIVE BROTH (CSB) (ISO 22964:2017)

### INSTRUCTION FOR USE

Rehydrate the contents of 1 vial with 5 ml of sterile distilled water and aseptically add 2 vials to 1000 ml of sterile molten TM 2047 - CRONOBACTER SELECTIVE BROTH (CSB) (ISO 22964:2017) / TM 1642 - VRE AGAR BASE (VANCOMYCIN RESISTANT ENTEROCOCCI AGAR) along with 1 vial of Meropenem Supplement to TM 1642. Mix well and pour into sterile Petri plates.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

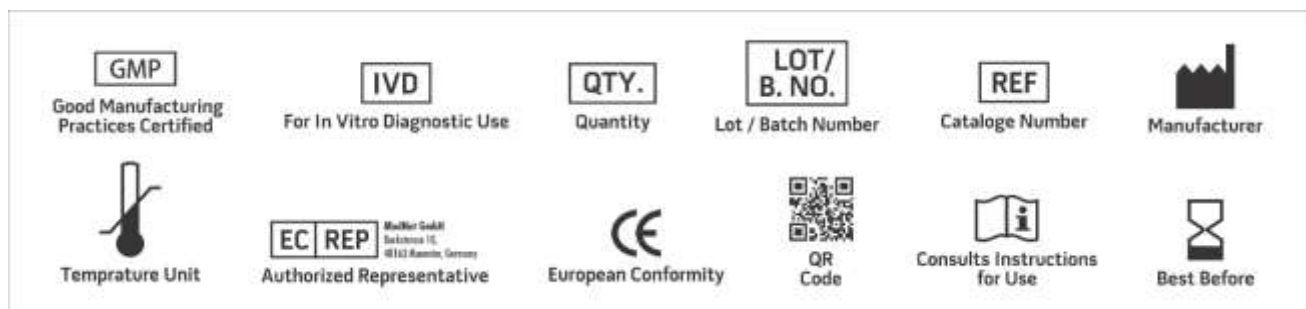

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

\*For Lab Use Only  
Revision: 2 March., 2022