

TS 266 – VITAMINS K1 SUPPLEMENT

INTENDED USE

A vitamin growth supplement used for the isolation of anaerobic organisms.

COMPOSITION

Ingredients	Concentration
Vitamin K1	10mg

(per vial sufficient for 1000 ml medium)

Applied in Medium


TM 1960 - ANAEROBIC BLOOD AGAR BASE

INSTRUCTION FOR USE

Rehydrate the content of 1 vial aseptically with 5 ml sterile distilled water. Mix well and aseptically add it to 1000 ml sterile, molten, cooled (45-50°C) TM 1960 - ANAEROBIC BLOOD AGAR BASE. Add 50 ml sterile defibrinated sheep blood to all the above media. Mix gently and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

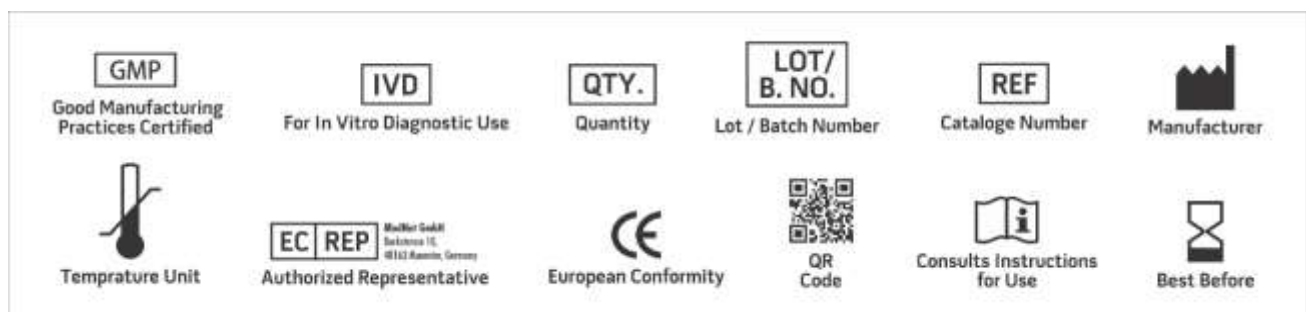

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 3 March., 2022