

TS 269 – BROMO THYMOL BLUE SUPPLEMENT

INTENDED USE

For the isolation and presumptive identification of microflora in urinary tract.

COMPOSITION

Ingredients	Concentration
Bromo thymol blue	20mg

(Per vial sufficient for 1000 ml medium)

Applied in Medium


TM 2025 - C.L.E.D. AGAR BASE W/O INDICATOR

INSTRUCTION FOR USE

Rehydrate the content of 1 vial aseptically with 2 ml of sterile 0.2 N Sodium hydroxide. Mix well and aseptically add to 1000 ml of sterile, molten, cooled (45-50°C) TM 2025 - C.L.E.D. AGAR BASE W/O INDICATOR. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

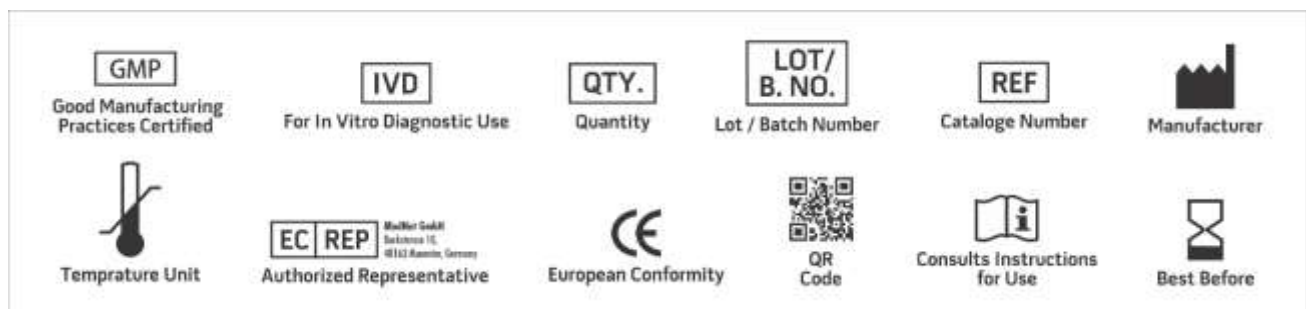

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 5 March., 2022