

## TS 300 – GROWTH SUPPLEMENT II FOR MSM

### INTENDED USE

A growth supplement for the enrichment of *Salmonella* species.

### COMPOSITION

| Ingredients  | Concentration |
|--------------|---------------|
| Ethanolamine | 2ml           |

(per vial sufficient for 500 ml medium)

### Applied in Medium


TM 2218 - MSM BROTH BASE

### INSTRUCTION FOR USE

Warm up the refrigerated contents of one vial to room temperature. Aseptically add the contents of one vial of Growth Supplement II for MSM (TS 300) to 1000 ml sterile, cooled (45-50°C) TM 2218 - MSM BROTH BASE along with Growth Supplement I for MSM (TM 299) . Mix well and dispense as desired.

### STORAGE

Vials should be stored in sealed container at 2-8°C.

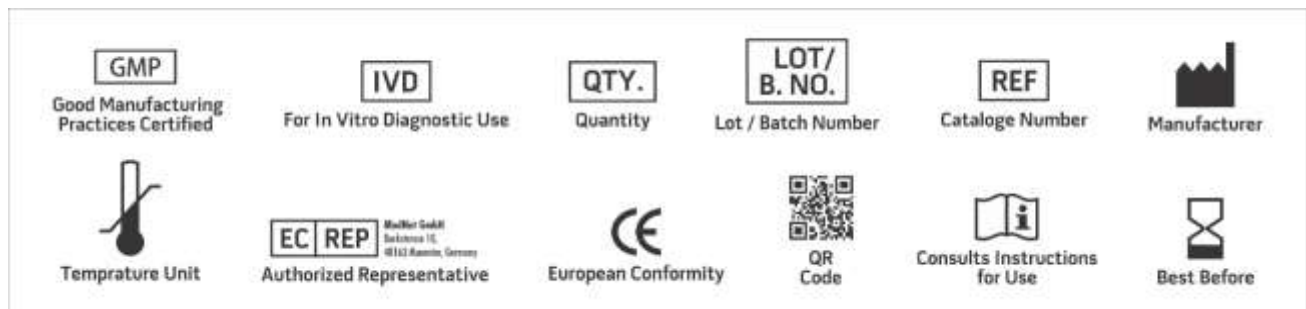

**NOTE:** Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

**\*For Lab Use Only**  
Revision: 28 Feb., 2022