

TS 320 – M-CP SELECTIVE SUPPLEMENT II

INTENDED USE

Recommended for selective isolation of *Clostridium perfringens*.

COMPOSITION

Ingredients	Concentration
Phenolphthalein diphosphate	0.050g
Distilled water	10ml

(Per vial sufficient for 500 ml medium)

Applied in Medium


TM 2130 - M-CP AGAR BASE

INSTRUCTION FOR USE

Warm up the refrigerated 0.5% phenolphthalein diphosphate solution to room temperature and add aseptically 10 ml of solution to 485 ml sterile, molten, cooled (45-50°C) TM 2130 - M-CP AGAR BASE along with rehydrated contents of one vial of M-CP Selective Supplement I / M-CP Selective Supplement, Modified. Mix well and pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

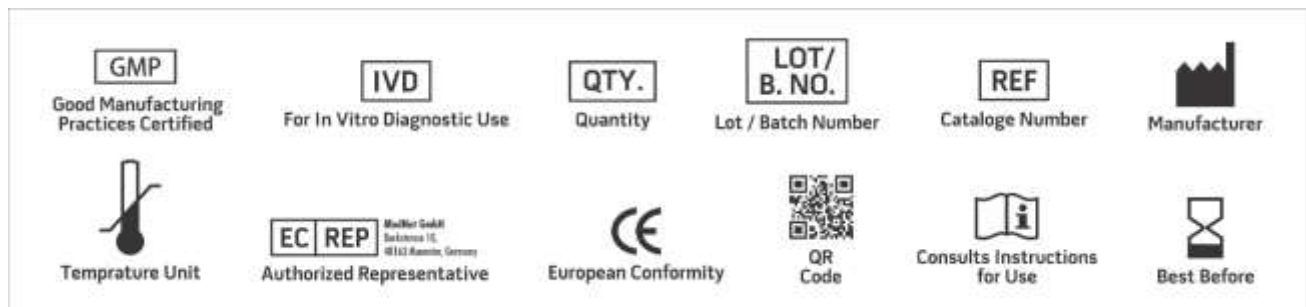

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 7 March., 2022